



ThermaSolutions

the original

Hyperthermic Perfusion Company

**ThermoChem™
HT-2000
&
HT-2500**

*Treating patients
worldwide since
our FDA approval
in 1999*

*Fully focused
and dedicated
to this
treatment*

*Easy to use
equipment and
disposables*



*Training and
Support play a
key role in our
strategy*

www.thermasolutions.com


The ThermoChem™ HT family of devices is the successor of the first fully-integrated system specifically designed and manufactured for Hyperthermic Perfusion. By exploiting heat and harnessing its power, this advanced medical technology offers a new choice in adjunctive surgical therapies. Used intraoperatively, the ThermoChem HT-2000 Series of devices raise the temperature to a pre-determined target, selected by the physician, of the peritoneal cavity by continuous warming and circulating of a sterile solution.

The ThermoChem™ HT System is optimized for an easy assembly of the disposable, quick heating, accurate and efficient use, safe waste management

The ThermoChem™ HT - 2000 Series devices includes two major components:

- **The ThermoChem™ HT-2000 & HT-2500* devices**
 - Adjustable temperature from 36 to 47°C in 0.1 increments
 - Flowrate of 0 to 2400 ml/min with adjustable roller pump
 - Consistent outflow with gravity and/or vacuum
 - Interactive touchscreen monitor for accurate readout or setting of patient fluid temperature, flowrate and system pressure levels
 - Data recording/downloading of patient procedure

- **The single-use Hyperthermic Perfusion Kit**
 - Optimized for both open and closed procedures
 - Packaged as an all-inclusive, pre-connected complete set with optional outflow tubes
 - 4L fluid reservoir
 - High efficient heat exchanger
 - Patient temperature probes, pressure sensor
 - Outflow tubing options:
 - Individual suction tubes or Lasso tubing for closed technique or Laparoscopic Lasso for HIPEC performed Laparoscopically
 - Inflow and outflow, multi-perforated silicon drains with outflow protectors
 - Incorporated waste bag for a completely closed system

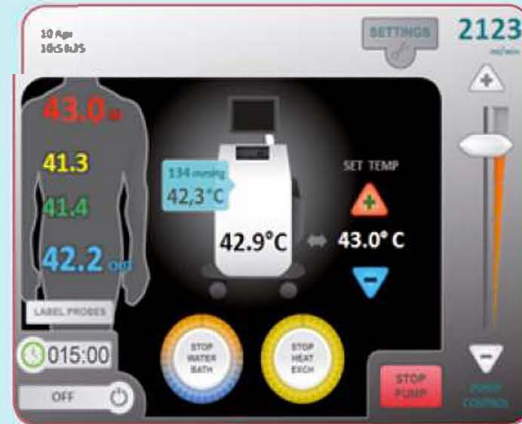

ThermaSolutions services the 5 continents from offices in both Europe and the USA:

ThermaSolutions LLC

1889 Buerkle Road
 White Bear Lake, MN 55110
 USA
 Tel. : +1 651 209 3900
 Fax : +1 651 209 3903
 E-mail : info@thermasolutions.com

ThermaSolutions Europe BV

Takkebijsters 41
 4817 BL Breda
 The Netherlands
 Tel. : +31 76 579 1144
 Fax : +31 76 579 1149
 E-mail : emea@thermasolutions.com

Support Program: Clinical support through reference centers, User / Service training and OR support.

*Check availability in your country.